

KAPA HTP Library Preparation Kit Illumina® platforms

KR0426 – v3.13

This Technical Data Sheet provides product information and a detailed protocol for the KAPA HTP Library Preparation Kit (Illumina® platforms), product codes KK8234 and KK8235.

Contents

Product Description	2
Product Applications	2
Product Specifications	2
Shipping and storage	2
Handling	2
Quality control	2
Important Parameters	3
Automated library construction	3
Reaction setup	3
Safe stopping points	3
Paramagnetic SPRI® beads and reaction cleanups ..	4
Input DNA and fragmentation	4
Cleanups after end repair and A-tailing	4
Adapter design and concentration	5
Post-ligation cleanup	5
Size selection	6
Library amplification	6
Evaluating the success of library construction	7
Process Workflow	9
Library Construction Protocol	10
Appendix 1:	
Library construction guidelines for target enrichment with the Roche Nimblegen™ SeqCap EZ system	17

Kit Codes and Components

KK8234 KK8235 96 libraries	KAPA End Repair Buffer (10X)	1.2 ml
	KAPA End Repair Enzyme Mix	600 µl
	KAPA A-Tailing Buffer (10X)	650 µl
	KAPA A-Tailing Enzyme	360 µl
	KAPA Ligation Buffer (5X)	1.3 ml
	KAPA DNA Ligase	600 µl
	KAPA PEG/NaCl SPRI® Solution	40 ml
	KAPA HiFi HotStart ReadyMix (2X)*	3 ml

*KK8235 is available for PCR-free workflows and does not contain KAPA HiFi HotStart ReadyMix for library amplification.

Quick Notes

- The protocol provided in this document is a generic prototype, and may require additional tailoring and optimization.
- The process workflow (p. 9) provides an overview of the library construction process and options for size selection.
- The **KAPA NGS Library Preparation Technical Guide** contains more detailed information about library construction parameters, and may facilitate protocol development and optimization.
- Separate, concentrated enzyme formulations and reaction buffers for end repair, A-tailing, and ligation provide the best combination of product stability, convenience, and efficiency.
- Adapters and PCR primers are not supplied with this kit, and can be obtained from any reputable oligonucleotide vendor.
- SPRI® beads are not included in the kit, but the PEG/NaCl SPRI® Solution required for "with-bead" reaction cleanups is provided.
- Generous reagent excesses are supplied to accommodate the dead volumes required for automated liquid handling.
- A single kit per batch of 96 samples simplifies reagent handling and inventory control.

Product Description

The KAPA HTP Library Preparation Kit is designed for high-throughput library construction for Illumina® sequencing, starting from fragmented, double-stranded DNA. The kit provides all of the enzymes and reaction buffers required for the following steps of library construction:

1. **End repair**, which produces blunt-ended, 5'-phosphorylated fragments.
2. **A-tailing**, during which dAMP is added to the 3'-ends of blunt-ended dsDNA library fragments.
3. **Adapter ligation**, during which dsDNA adapters with 3'-dTTP overhangs are ligated to 3'-A-tailed library fragments.
4. **Library amplification (optional)**, which employs PCR to amplify library fragments carrying appropriate adapter sequences on both ends.

The kit has been validated for library construction from 100 ng – 5 µg of human genomic DNA for whole-genome shotgun sequencing or targeted sequencing by solution hybrid selection (capture). For smaller genomes, or lower complexity samples, such as ChIP DNA, amplicons, or cDNA (for RNA-seq), successful library construction has been achieved from low nanogram to picogram quantities (≥100 pg) of input DNA.

The kit provides all of the enzymes and buffers required for library construction and amplification, but does not include adapters, PCR primers or SPRI® beads. Enzyme formulations and reaction buffers for end repair, A-tailing and ligation are supplied in convenient, concentrated formats.

Efficient, cost-effective and automation-friendly reaction cleanups and higher recovery of input DNA are achieved through implementation of the "with-bead" strategy developed at The Broad Institute of MIT & Harvard and Foundation Medicine¹. The kit includes PEG/NaCl SPRI® (Solid Phase Reversible Immobilization) Solution for this purpose.

In order to maximize sequence coverage uniformity, it is critical to minimize library amplification bias. KAPA HiFi DNA Polymerase has been designed for low-bias, high-fidelity PCR, and is the reagent of choice for NGS library amplification^{2, 3, 4}. The KAPA HTP Library Preparation Kit (KK8234) includes KAPA HiFi HotStart ReadyMix (2X), a ready-to-use PCR mix comprising all the components for library amplification, except primers and template. A kit without an amplification module (KK8235) is available for PCR-free workflows. These kits can also be combined with KAPA HiFi Real-Time Library Amplification Kits (KK2701 and KK2702), or with KAPA HiFi HotStart Uracil+ ReadyMix (KK2801 and KK2802) for the amplification of libraries that have undergone bisulfite-treatment.

1. Fisher, S. *et al. Genome Biology* **12**, R1 (2011).
2. Oyola, S.O. *et al. BMC Genomics* **13**, 1 (2012).
3. Quail M.A. *et al. Nature Methods* **9**, 10 – 11 (2012).
4. Quail M.A. *et al. BMC Genomics* **13**: 341 (2012).

Product Applications

The KAPA HTP Library Preparation Kit is ideally suited for high-throughput NGS library construction workflows that involve end repair, A-tailing, adapter ligation, and library amplification (optional). The kit has been tailored to accommodate the specific requirements of automated liquid handling, and the protocol may be adapted for incorporation into workflows for a wide range of NGS applications, including:

- Whole-genome shotgun sequencing
- Targeted sequencing by solution hybrid selection (i.e. exome or custom capture using the Roche Nimblegen™, Agilent SureSelect, Illumina® TruSeq™, or IDT xGen™ Lockdown™ Probes systems)
- ChIP-seq
- RNA-seq
- Methyl-seq (in combination with the KAPA HiFi HotStart Uracil+ ReadyMix) Library Amplification ReadyMix)

Specific guidelines for the construction of libraries for target enrichment using the Roche Nimblegen™ SeqCap EZ system may be found in Appendix 1.

Product Specifications

Shipping and storage

The enzymes provided in this kit are temperature sensitive, and appropriate care should be taken during shipping and storage. KAPA Library Preparation Kits are shipped on dry ice or ice packs, depending on the country of destination. Upon receipt, immediately store enzymes and reaction buffers (including PEG/NaCl SPRI® Solution) at -20 °C in a constant-temperature freezer. When stored under these conditions and handled correctly, the kit components will retain full activity until the expiry date indicated on the kit label.

Handling

Always ensure that components have been fully thawed and thoroughly mixed before use. Keep all enzyme components and master mixes on ice as far as possible during handling and preparation. KAPA HiFi HotStart ReadyMix (2X) contains isostabilizers and may not freeze completely, even when stored at -20 °C. Nevertheless, always ensure that the KAPA HiFi HotStart ReadyMix is fully thawed and thoroughly mixed before use. PEG/NaCl SPRI® Solution does not freeze at -20 °C, but should be equilibrated to room temperature and thoroughly mixed before use.

Quality control

All kit components are subjected to stringent functional quality control, are free of detectable contaminating exo- and endonuclease activities, and meet strict requirements with respect to DNA contamination. Please contact support@kapabiosystems.com for more information.

Important Parameters

High-throughput library construction workflows must be tailored and optimized to accommodate specific experimental designs, sample characteristics, sequencing applications, and equipment. The protocol provided in this document is a generic prototype, and there are many parameters which may be adjusted to optimize performance, efficiency, and cost-effectiveness.

In addition to the information in this section, please consult the **KAPA NGS Library Preparation Technical Guide** and/or contact support@kapabiosystems.com for further guidelines when designing or optimizing your library construction workflow.

Automated library construction

The library construction protocol described in this document can be carried out manually, and most protocol development and validation work is usually done this way. Although it is possible to achieve moderately high sample throughput by using multi-channel pipettes and 96-well plates, automated liquid handling is indispensable for high-throughput NGS production lines, and automating a validated manual library construction protocol can represent a significant challenge.

In addition to increased sample throughput, automation may be expected to provide additional advantages such as improved reproducibility and process control. Nevertheless, automation may result in slightly lower yields and/or different size distributions when compared with manual library construction performed by a skilled, experienced and attentive technician. Most often, these discrepancies can be minimized through careful selection of appropriate hardware and plasticware, and by optimizing liquid handling parameters such as aspiration speeds and volumes in automation scripts.

Kapa Biosystems does not supply automated liquid handling equipment, but we are constantly working in partnership with automation solution providers and customers to develop and validate optimized, automated methods for liquid handling platforms suitable for use in NGS library construction. Please contact support@kapabiosystems.com for more information about using this kit with your particular automated liquid handling system.

Reaction setup

This kit is intended for high-throughput library construction, and the protocol is therefore designed to be automation-friendly. For this reason, and to enable a streamlined “with-bead” strategy, reaction components are combined into master mixes, rather than dispensed separately into individual reactions. When processing multiple samples, prepare 10 – 20% excess of each master mix. When processing batches of 48 or more samples with an automated liquid handling system, 20% excess is required for all master mixes. For other

reagents (adapters, AMPure® XP reagent, SPRI® Solution, 80% ethanol and elution buffer), the required excess may vary from one specific liquid handling system to another. Please refer to Section 1 of the Protocol for more details.

96-well PCR plates are typically used for high-throughput library construction. The maximum working volume in these plates is usually ~200 µl, and this is accounted for in this protocol. It may be possible to employ 96-well plates with larger working volumes or deep well plates to accommodate larger reaction volumes for special applications. Always use plastics that are certified to be free of DNAses, RNAses, and nucleases. Low DNA-binding plastics are recommended. When selecting the most appropriate plasticware for your workflow, consider compatibility with:

- the plate gripper and other components of your liquid handling system.
- the magnet used during SPRI® bead manipulations.
- vortex mixers and centrifuges, where appropriate.
- heating blocks or thermocyclers used for reaction incubations and/or library amplification.

Safe stopping points

The library construction process, from end repair to final, amplified library, can be performed in 4 – 8 hours, depending on the specific workflow and number of samples being processed. Automated methods are typically designed to complete the process from end repair to post-ligation processing in an uninterrupted manner, with minimal user intervention. However, the protocol may be paused safely after any of the bead cleanup steps, as described below:

- After the end repair cleanup (Steps 3.1 – 3.13), resuspend the washed beads in 20 µl of 1X A-Tailing Buffer (without enzyme; Step 4.1), and store the reactions at 4 °C for up to 24 hours.
- After the A-tailing cleanup (Steps 5.1 – 5.13), resuspend the washed beads in 20 µl of 1X Ligation Buffer (without enzyme or adapter; Step 6.1), and store the reactions at 4 °C for up to 24 hours.
- After the first post-ligation cleanup (Steps 7.1 – 7.13), resuspend the washed beads in the appropriate volume of 10 mM Tris-HCl (pH 8.0) as outlined in Step 7.14, and store the reactions at 4 °C for up to 24 hours.

DNA solutions containing beads must not be frozen, and beads must not be stored dry, as this is likely to damage the beads and result in sample loss. To resume the library construction process, centrifuge the reaction vessels briefly to recover any condensate, and add the remaining components required for the next enzymatic reaction in the protocol (see Tables 4B and 5B on p. 10). If the protocol was paused after the first post-ligation cleanup, continue directly with the second post-ligation cleanup (Step 7.16), dual-SPRI® size selection (Step 8.1), or size selection using an alternative method.

Safe stopping points (continued)

Adapter-ligated DNA that has been completely cleaned up or size-selected can be stored at 4 °C for one week, or at -20 °C for at least one month before amplification, target enrichment, and/or sequencing. To avoid degradation, always store DNA in a buffered solution (10 mM Tris-HCl, pH 8.0) and minimize the number of freeze-thaw cycles.

Paramagnetic SPRI® beads and reaction cleanups

- Cleanups should be performed in a timely manner to ensure that enzyme reactions do not proceed beyond optimal incubation times.
- This protocol has been validated using Agencourt® AMPure® XP reagent (Beckman Coulter, part number A63880, A63881, or A63882). Solutions and conditions for DNA binding and size selection may differ if other beads are used.
- Observe all the manufacturer's storage and handling recommendations for AMPure® XP reagent.
- Beads will settle gradually; always ensure that they are fully resuspended before aspirating AMPure® XP reagent.
- The incubation times provided for reaction cleanups and size selection are guidelines only, and should be modified/optimized according to your current protocols, previous experience, and specific equipment and samples in order to maximize library construction efficiency and throughput.
- The time required to completely capture magnetic beads varies according to the reaction vessel and magnet used. It is important to not discard or transfer any beads with the removal or transfer of supernatant. Capture times should be optimized accordingly.
- The volumes of 80% ethanol used for bead washes may be adjusted to accommodate smaller reaction vessels and/or limited pipetting capacity, but it is important that the beads are entirely submerged during the wash steps. Where possible, use a wash volume that is equal to the volume of sample plus AMPure® XP reagent or PEG/NaCl SPRI® Solution.
- It is important to remove all ethanol before proceeding with subsequent reactions. However, over-drying of beads may make them difficult to resuspend, and may result in a dramatic loss of DNA. With optimized pipetting, drying of beads for 3 – 5 min at room temperature should be sufficient. **Drying of beads at 37 °C is not recommended.**
- Where appropriate, DNA should be eluted from beads in elution buffer (10 mM Tris-HCl, pH 8.0). Elution of DNA in PCR-grade water is not recommended, as DNA is unstable in unbuffered solutions.

Input DNA and fragmentation

- This protocol has been validated for library construction from 100 ng – 5 µg of appropriately fragmented, double-stranded DNA. However, libraries can be prepared from lower input amounts if the sample represents sufficient copies to ensure the requisite coverage and complexity in the final library. Successful library construction has been achieved from <100 pg of ChIP DNA, low nanogram quantities of cDNA or microbial DNA, and 1 – 10 ng of high-quality human or mouse genomic DNA.
- The above typically refers to the input into the end repair reaction. If input DNA is quantified before fragmentation, and/or fragmented DNA is subjected to cleanup or size selection prior to end repair, the actual input into library construction may be significantly lower. This should be taken into account when evaluating the efficiency of the process and/or during optimization of library amplification cycle number.
- The proportion of fragmented DNA that is successfully converted to adapter-ligated molecules decreases as input is reduced. When starting library construction (end repair) with >100 ng fragmented DNA, 15 – 40% of input DNA is typically recovered as adapter-ligated molecules, whereas the recovery typically ranges from 0.5 to 15% for libraries constructed from 100 pg – 100 ng DNA. These figures apply to high quality DNA and can be significantly lower for DNA of lower quality, e.g. FFPE samples. Workflows that contain additional SPRI® cleanups or size selection prior to library amplification are likely to result in a lower yield of adapter-ligated molecules.
- Solutions containing high concentrations of EDTA and strong buffers may negatively affect the end repair reaction, and should be avoided. If fragmented DNA will not be processed (i.e. subjected to cleanup or size selection) prior to end repair, DNA should be fragmented in 10 mM Tris-HCl (pH 8.0 or 8.5) with 0.1 mM EDTA. Fragmentation in water is not recommended.
- In some circumstances it may be convenient to fragment input DNA in 1X KAPA End Repair Buffer, in which case the end repair reaction setup should be adjusted accordingly. Please contact support@kapabiosystems.com for more information.

Cleanups after end repair and A-tailing

- This protocol provides for 1.7X – 1.8X cleanups after end repair and A-tailing. This ratio of PEG/NaCl SPRI® Solution to sample volume will retain most DNA fragments larger than ~75 bp. If you wish to retain very small DNA fragments, the PEG/NaCl SPRI® Solution to sample ratio can be increased to 2X – 3X for all cleanups prior to adapter ligation.

Cleanups after end repair and A-tailing (continued)

- If a >2X SPRI® bead cleanup is desired after end repair, the end repair reaction must be scaled down when performing library construction in standard PCR plates (maximum working volume of ~200 ul). Please contact support@kapabiosystems.com if your workflow or sample type requires modified SPRI® bead cleanups.

Adapter design and concentration

- This protocol has been validated using standard, indexed illumina® TruSeq™ "forked" adapters, but the kit is compatible with other adapters of similar design.
- Adapter concentration affects ligation efficiency, as well as adapter and adapter-dimer carry-over in post-ligation cleanups. The optimal adapter concentration for your workflow represents a compromise between cost and the above factors. Your choice of post-ligation cleanup and size-selection options should be informed by your choice of adapter concentration. Please refer to **Important Parameters: Post-ligation cleanup** for more details.
- Ligation efficiency is robust for adapter:insert molar ratios ranging from 10:1 to 50:1. As a general guideline, we recommend an adapter:insert molar ratio of ~10:1, for libraries constructed from ≥100 ng fragmented DNA. This translates to different final adapter concentrations for libraries with different size distributions (see Table 1 below). An adapter:insert molar ratio >10:1 may be beneficial for libraries constructed from lower amounts of input DNA.
- While it is not necessary to adjust adapter concentrations to accommodate moderate sample-to-sample variations, we recommend using an adapter concentration that is appropriate for the range of input DNA concentrations.
- The best way to accommodate different adapter concentrations within a batch of samples processed together, is to vary the concentration of adapter stock solutions, and dispense a fixed volume (5 µl) of each

adapter. The alternative – using a single stock solution, and dispensing variable volumes of adapter into ligation reactions – is not recommended for automated workflows.

Post-ligation cleanup

- It is important to remove excess unligated adapter and adapter-dimer molecules from the library prior to library amplification or cluster generation.
- While a single SPRI® bead cleanup removes most unligated adapter and adapter-dimer, a second SPRI® bead cleanup is recommended to eliminate any remaining adapter species from the library. The amount of adapter and adapter-dimer carried through the first cleanup is dependent on the adapter concentration in the ligation reaction.
- If size selection is carried out between adapter ligation and library amplification (or clustering), a single post-ligation cleanup with SPRI® beads (1X) is usually sufficient prior to size selection. If no post-ligation size selection is carried out, two consecutive 1X SPRI® bead cleanups are recommended.
- The volume in which washed beads are resuspended after the post-ligation cleanup(s) should be adjusted to suit your chosen workflow:
 - If proceeding directly to library amplification, determine an appropriate final volume in which to elute the library DNA, keeping in mind that you may wish to divert and/or reserve some of this library material for archiving and/or QC purposes. Since an optimized 50 µl library amplification reaction should yield ~1 µg of DNA, and can accommodate a maximum of 20 µl template DNA, an elution volume of 22 – 32 µl is recommended.
 - If proceeding with size selection, elute the library DNA in an appropriate volume according to the size selection method of choice. For the dual-SPRI® size selection procedure described here, beads have to be resuspended in a final elution volume of 100 µl.

Table 1. Recommended adapter concentrations.

Insert DNA per 50 µl end repair reaction	Recommended adapter concentration for DNA sheared to an average size of					
	175 bp		350 bp		500 bp	
	Stock	Final	Stock	Final	Stock	Final
3 – 5 µg	60 µM	6 µM	30 µM	3 µM	21 µM	2.1 µM
1 µg	20 µM	2 µM	10 µM	1 µM	7 µM	0.7 µM
500 ng	10 µM	1 µM	5 µM	500 nM	3.5 µM	350 nM
100 ng	2 µM	200 nM	1 µM	100 nM	700 nM	70 nM
10 ng	200 nM	20 nM	100 nM	10 nM	70 nM	7 nM

Size selection

- Size selection requirements vary widely according to specific applications. Depending on preference, the dual-SPRI® size selection procedures presented in this protocol may be omitted entirely, modified, or replaced with alternative size selection procedures. Size selection may be carried out at alternative points in the overall workflow, for example:
 - prior to end repair of fragmented DNA.
 - immediately before library amplification (as outlined in this protocol).
 - after library amplification.
- Size selection inevitably leads to a loss of sample material. Depending on the details, these losses can be dramatic (>80%), and may significantly increase the number of amplification cycles needed to generate sufficient material for the next step in the process (capture or sequencing). The potential advantages of one or more size selection steps in a library construction workflow should be weighed against the potential loss of library complexity, especially when input DNA is limited. A carefully optimized fragmentation protocol, especially for shorter insert libraries and/or read lengths, may eliminate the need for size selection, thereby simplifying the library construction process and limiting sample losses.
- Over-amplification of libraries often results in the observation of secondary, higher molecular weight peaks when amplified libraries are analyzed electrophoretically. These higher molecular weight peaks are artefacts, and typically contain authentic library molecules of the appropriate length. To eliminate these artefacts, optimization of library amplification reaction parameters (cycle number and/or primer concentration), rather than post-amplification size selection, is recommended. Please refer to the **Library amplification** subsection for more information.
- KAPA Ligation Buffer contains high concentrations of PEG 6000, which will interfere with efficient dual-SPRI® size selection and can affect the efficiency of other size selection techniques if not removed. It is therefore important to perform at least one post-ligation SPRI® bead cleanup (1X) prior to dual-SPRI® or any other size selection method.
- The dual-SPRI® size selection procedure described in Section 8 of the protocol is designed for selection of adapter-ligated fragments approximately 250 – 450 bp in size. Consult the **KAPA NGS Library Preparation Technical Guide** if you wish to select a different range of fragment sizes.
- "Forked" adapters with long single-stranded arms (such as the TruSeq™ design) noticeably affect the size-dependent binding of DNA to SPRI® beads, as well as the apparent size of fragments determined by some electrophoresis instruments (e.g. those

employing microfluidic chips and tapes). Size selection parameters will therefore require optimization depending on a number of factors, including:

- The design of the adapters used.
- The point at which size selection is applied in the protocol; size selection of dsDNA fragments prior to library construction or after amplification may require different parameters than those used for post-ligation size selection of fragments carrying forked adapter ends.
- Dual-SPRI® size selection is sensitive to multiple factors that are beyond the scope of this document. The **KAPA NGS Library Preparation Technical Guide** contains additional guidelines for the optimization of dual-SPRI® size selection parameters. Any dual-SPRI® size selection protocol should be carefully optimized and validated before it is used for precious samples.

Library amplification

- KAPA HiFi HotStart, the enzyme provided in the KAPA HiFi HotStart ReadyMix, is an antibody-based hot start formulation of KAPA HiFi DNA Polymerase, a novel B-family DNA polymerase engineered for increased processivity and high fidelity. KAPA HiFi HotStart DNA Polymerase has 5'→3' polymerase and 3'→5' exonuclease (proofreading) activities, but no 5'→3' exonuclease activity. The strong 3'→5' exonuclease activity results in superior accuracy during DNA amplification. The error rate of KAPA HiFi HotStart DNA Polymerase is 2.8×10^{-7} errors/base, equivalent to 1 error in 3.5×10^6 nucleotides incorporated.
- In library amplification reactions (set up according to the recommended protocol), primers are typically depleted before dNTPs. When DNA synthesis can no longer take place due to substrate depletion, subsequent rounds of DNA denaturation and annealing result in the separation of complementary DNA strands, followed by imperfect annealing to non-complementary partners. This presumably results in the formation of so-called "daisy-chains" or "tangled knots", comprising large assemblies of improperly annealed, partially double-stranded, heteroduplex DNA. These species migrate slower and are observed as secondary, higher molecular weight peaks during the electrophoretic analysis of amplified libraries. However, they are typically comprised of library molecules of the desired length, which are individualized during denaturation prior to cluster amplification or probe hybridization. Since these heteroduplexes contain significant portions of single-stranded DNA, over-amplification leads to the under-quantification of library molecules with assays employing dsDNA-binding dyes. qPCR-based library quantifications methods, such as the KAPA Library Quantification assay, quantify DNA by denaturation and amplification, thereby providing an

Library amplification (continued)

accurate measure of the amount of adapter-ligated molecules in a library, even if the library was over-amplified.

Please refer to the **KAPA NGS Library Preparation Technical Guide** for a more detailed discussion of factors that can affect the efficiency of library amplification, and the impact of over-amplification on library quantification.

- Excessive library amplification can result in other unwanted artefacts such as amplification bias, PCR duplicates, chimeric library inserts, and nucleotide substitutions. The extent of library amplification should therefore be limited as much as possible, while ensuring that sufficient material is generated for QC and downstream processing (e.g. target enrichment or sequencing).
- If cycled to completion (**not recommended**) a single 50 µl library amplification PCR, performed as described in Section 9, can produce 8 – 10 µg (160 – 200 ng/µl) of amplified library. To minimize over-amplification and associated undesired artefacts, the number of amplification cycles should be optimized to produce an amplified library with a concentration in the range of 10 – 30 ng/µl (0.5 – 1.5 µg of PCR product per 50 µl reaction).
- Quantification of adapter-ligated libraries prior to library amplification can greatly facilitate the optimization of library amplification parameters, particularly when a library construction workflow is first established. With the KAPA Library Quantification Kit, the amount of template DNA (adapter-ligated molecules) available for library amplification can be determined accurately. Using a simple calculation

Table 2. Theoretical number of cycles required to achieve 0.5 – 1.5 µg (10 – 30 ng/µl) library in a standard 50 µl KAPA HiFi HotStart ReadyMix library amplification reaction, starting from different amounts of template DNA.

Template DNA	Number of cycles
1 ng	10 – 11
5 ng	7 – 8
10 ng	6 – 7
25 ng	5 – 6
50 ng	4 – 5
100 ng	3 – 4
250 ng	2 – 3
500 ng	0 – 2

(for exponential amplification), the number of amplification cycles needed to achieve a specific yield of amplified library may be predicted theoretically (Table 2).

- The **actual** optimal number of amplification cycles may be 1 – 3 cycles higher, particularly for libraries constructed from FFPE DNA or other challenging samples, or libraries with a broad fragment size distribution.


Evaluating the success of library construction

- Your specific library construction workflow should be tailored and optimized to yield a sufficient amount of adapter-ligated molecules of the desired size distribution for the next step in the process (e.g. target enrichment or sequencing), as well as for QC and archiving purposes.
- The size distribution of the final or pre-capture library should be confirmed with an electrophoretic method, whereas KAPA Library Quantification Kits for Illumina® platforms are recommended for qPCR-based quantification of libraries. These kits employ primers based on the Illumina® flow cell oligos, and can be used to quantify libraries:
 - that are ready for flow-cell amplification.
 - that were constructed with full-length adapters after ligation, size selection, and pre-capture amplification.
- Once a library construction workflow has been optimized, and consistently yields the desired amount of amplified library, it is typically not necessary to perform in-process quality control. However, qPCR-based quantification of libraries after adapter ligation or prior to library amplification can provide useful data for optimization or troubleshooting.
- Accurate quantification of DNA at this stage will allow you to evaluate the efficiency of:
 - the process from end repair to ligation, by determining the percentage of input DNA converted to adapter-ligated molecules.
 - library amplification with the selected number of cycles, based on the actual amount of template DNA used in the PCR.
- The availability of quantification data before and after library amplification allows the two major phases of the library construction process to be evaluated and optimized independently to achieve the desired yield of amplified library.
- If size selection is performed at any stage of the process, qPCR quantification before and after size selection may also be helpful to define the relative benefit of size selection, and the loss of material associated with the process.

Evaluating the success of library construction (continued)

- Electrophoretic evaluation of libraries after adapter ligation or before library amplification may be informative, but remember that unamplified libraries prepared with "forked" adapters will appear to have a larger, wider or bimodal fragment size distribution than the corresponding amplified libraries (refer to **Size selection** on p. 6). The difference in overall appearance and fragment size distribution of an unamplified vs corresponding amplified library varies depending on the electrophoretic system used. To accurately determine the size distribution of an unamplified library, an aliquot of the library may be subjected to one or two cycles of amplification (to ensure that adapters are fully double-stranded) prior to electrophoretic analysis. Alternatively, useful information may also be obtained by electrophoretic analysis of library quantification products generated with the KAPA Library Quantification Kit.

Process Workflow

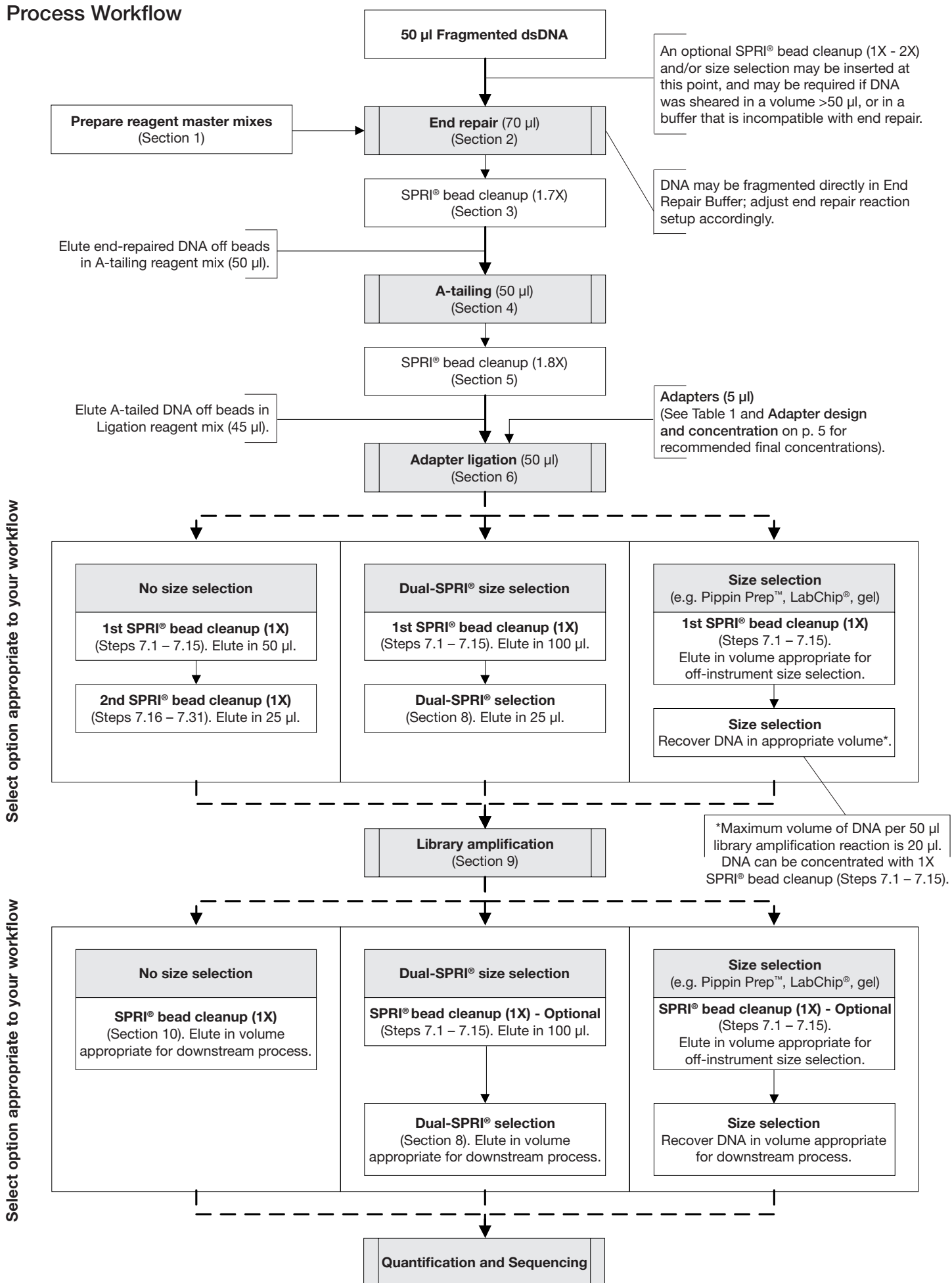

Library Construction Protocol

1. Reagent preparation

For maximum stability and shelf-life, the enzyme formulations and concentrated reaction buffers for end repair, A-tailing and ligation are supplied separately in KAPA HTP Library Preparation Kits. For a streamlined "with-bead" protocol, a reagent master mix is prepared for each of these enzymatic steps, as outlined in Tables 3 – 6.

Master mixes may be constituted with varying proportions of the total final water requirement. In the examples given in the tables below, all the required water is included in each master mix, allowing the entire reaction mix to be added in a single pipetting step (after cleanups, beads are resuspended directly in master mix for the next enzymatic step). This strategy may be modified for some automated liquid handling systems, i.e. master mixes may be prepared with only a portion of the total volume of water, and the rest of the water added directly to the beads, or to beads resuspended in the reaction mix.

At safe stopping points, some or all of the water and/or reaction buffer can be added to the beads, for storage at 4 °C for up to 24 hours. To resume library construction, prepare the master mix with the remaining volume of water and reaction buffer, and the required volume of enzyme. Recommendations on how to formulate master mixes at safe stopping points are provided in Tables 4B and 5B.

Master mixes for end repair, A-tailing and ligation may be prepared immediately before use, or stored for up to 1 week at 4 °C or 2 weeks at -20 °C. Master mixes stored at -20 °C are stable through three freeze-thaw cycles.

Additional reagents required for the KAPA HTP library preparation protocol are listed in Table 7.

The excess of each master mix or reagent required for processing a full batch of 96 samples is given in each of the tables. This is based on current experience with different liquid handling platforms, and varies as a result of the way reagents are utilized by different systems. When processing batches of less than 48 samples, the maximum excess may not be required.

Table 3. Composition of the end repair reaction

Component	1 Library	Excess required	96 Libraries (incl. max. % excess)
End Repair Master Mix:			
Water	8 µl	5 – 20%	922 µl
10X KAPA End Repair Buffer	7 µl		806 µl
KAPA End Repair Enzyme Mix	5 µl		576 µl
Total master mix volume	20 µl		2 304 µl
Other components of the final reaction:			
Fragmented DNA	50 µl		
Balance of water required	0 µl		
Total reaction volume	70 µl		

Table 4A. Composition of the A-tailing reaction (standard, uninterrupted protocol)

Component	1 Library	Excess required	96 Libraries (incl. max. % excess)
A-Tailing Master Mix:			
Water	42 µl	5 – 20%	4 838 µl
10X KAPA A-Tailing Buffer	5 µl		576 µl
KAPA A-Tailing Enzyme	3 µl		346 µl
Total master mix volume	50 µl		5 760 µl
Other components of the final reaction:			
End repaired DNA with beads	0 µl		
Balance of water required	0 µl		
Total reaction volume	50 µl		

Table 4B. Composition of the A-tailing reaction if process is stopped after end repair.

Component	1 Library	Excess required	96 Libraries (incl. max. % excess)
Resuspend beads in mix of: Water 10X KAPA A-Tailing Buffer KAPA A-Tailing Enzyme Total volume of mix	18 µl 2 µl 0 µl 20 µl	5 – 20%	2 074 µl 230 µl 0 µl 2 304 µl
To resume, add mix of: Water 10X KAPA A-Tailing Buffer KAPA A-Tailing Enzyme Total volume of mix	24 µl 3 µl 3 µl 30 µl		2 765 µl 346 µl 346 µl 3 456 µl
Other components of the final reaction: End repaired DNA with beads Balance of water required Total reaction volume	0 µl 0 µl 50 µl		

Table 5A. Composition of the ligation reaction (standard, uninterrupted protocol).

Component	1 Library	Excess required	96 Libraries (incl. max. % excess)
Ligation Master Mix: Water 5X KAPA Ligation Buffer KAPA T4 DNA Ligase Total master mix volume	30 µl 10 µl 5 µl 45 µl	5 – 20%	3 456 µl 1 152 µl 576 µl 5 184 µl
Other components of the final reaction: A-tailed DNA with beads Adapter (<0.7 – 60 µM, as appropriate) Balance of water required Total reaction volume	0 µl 5 µl 0 µl 50 µl	50 – 100%	7.5 – 10 µl adapter per library

Table 5B. Composition of the ligation reaction if process is stopped after A-tailing.

Component	1 Library	Excess required	96 Libraries (incl. max. % excess)
Resuspend beads in mix of: Water 5X KAPA Ligation Buffer KAPA T4 DNA Ligase Total volume of mix	16 µl 4 µl 0 µl 20 µl	5 – 20%	1 843 µl 461 µl 0 µl 2 304 µl
To resume, add mix of: Water 5X KAPA Ligation Buffer KAPA T4 DNA Ligase Total volume of mix	14 µl 6 µl 5 µl 25 µl		1 613 µl 691 µl 576 µl 2 880 µl
Other components of the final reaction: A-tailed DNA with beads Adapter (<0.7 – 60 µM, as appropriate) Balance of water required Total reaction volume	0 µl 5 µl 0 µl 50 µl	50 – 100%	7.5 – 10 µl adapter per library

Table 6. Composition of the library amplification reaction.

Component	1 Library	Excess required	96 Libraries (incl. max. % excess)
Library Amplification Master Mix: Water 2X KAPA HiFi HotStart ReadyMix Total master mix volume	0 µl 25 µl 25 µl	5 – 20%	0 µl 2 880 µl 2 880 µl
Other components of the final reaction: Adapter-ligated library DNA PCR Primer Premix (5 µM each primer) Balance of water required Total reaction volume	20 µl 5 µl 0 µl 50 µl	5 – 20%	576 µl

Table 7. Approximate volumes of additional reagents required.

Reagent	1 Library	Excess required	96 Libraries (incl. max. % excess)
PEG/NaCl SPRI® Solution (provided in kit): A-tailing cleanup 1 st Adapter ligation cleanup 2 nd Adapter ligation cleanup/Dual-SPRI® size selection (max.) Total volume required	90 µl 50 µl 60 µl 200 µl	5 – 30%	260 µl per library or 25 ml for 96 libraries
AMPure® XP reagent (not supplied): End repair cleanup Dual-SPRI® size selection Post amplification cleanup Total volume required	120 µl ~20 µl 50 µl ~190 µl	5 – 20%	~230 µl per library or ~22 ml for 96 libraries
80% Ethanol (freshly prepared; not supplied): End repair cleanup A-tailing cleanup 1 st Adapter ligation cleanup 2 nd Adapter ligation cleanup/Dual-SPRI® size selection Post amplification cleanup Total volume required*	400 µl 400 µl 400 µl 400 µl 400 µl 2 ml	5 – 40%	2.8 ml per library or 270 ml for 96 libraries
Elution buffer (10 mM Tris-HCl, pH 8.0; not supplied): Adapter ligation cleanup/Dual-SPRI® size selection (max.) Amplification cleanup Total volume required	~125 µl ~25 µl ~150 µl	10 – 30%	195 µl per library or 18.7 ml for 96 libraries

*When using automated liquid handling systems with filtered pipette tips, the volume of 80% ethanol used per wash is often reduced to 150 µl to accommodate the maximum working volume of the tip. In such cases, the total volume of 80% ethanol required for per sample is 1.5 ml, or 2.1 ml with the maximum excess, which equates to ~200 ml for 96 samples.

2. End Repair Reaction Setup

2.1 Assemble each end repair reaction as follows:

Fragmented, double-stranded DNA	50 µl
End Repair Master Mix (Table 3)	20 µl
Total reaction volume	70 µl

2.2 Mix, and incubate at 20 °C for 30 min. Proceed immediately to the next step.

3. End Repair Cleanup

3.1 To each 70 µl end repair reaction, add:

Agencourt® AMPure® XP reagent	120 µl
Total volume per well/tube	190 µl

3.2 Mix thoroughly by pipetting up and down multiple times.

3.3 Incubate the plate/tube at room temperature for 5 – 15 min to allow the DNA to bind to the beads.

3.4 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.

3.5 Carefully remove and discard the supernatant.

3.6 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.

3.7 Incubate the plate/tube at room temperature for ≥30 sec.

3.8 Carefully remove and discard the ethanol.

3.9 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.

3.10 Incubate the plate/tube at room temperature for ≥30 sec.

3.11 Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.

3.12 Allow the beads to dry at room temperature, sufficiently for all the ethanol to evaporate. **Caution: over-drying the beads may result in dramatic yield loss.**

3.13 Remove the plate/tube from the magnet.

4. A-Tailing Reaction Setup

4.1 To each well/tube containing beads with end-repaired DNA, add:

Water (if required)	0 µl
A-Tailing Master Mix (Table 4A) [†]	50 µl
Total reaction volume	50 µl

[†]For a safe stopping point, resuspend the beads in 20 µl 1X A-Tailing Buffer, cover the reaction and store at 4 °C for up to 24 hours. Resume the reaction by adding the rest of the reaction components, as outlined in Table 4B (p. 11).

4.2 Thoroughly resuspend the beads by pipetting up and down multiple times.

4.3 Incubate at 30 °C for 30 min. Proceed immediately to the next step.

5. A-Tailing Cleanup

5.1 To each well/tube containing the 50 µl A-tailing reaction with beads, add:

PEG/NaCl SPRI® Solution	90 µl
Total volume per well/tube	140 µl

5.2 Thoroughly resuspend the beads by pipetting up and down multiple times.

5.3 Incubate the plate/tube at room temperature for 5 – 15 min to allow the DNA to bind to the beads.

5.4 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.

5.5 Carefully remove and discard the supernatant.

5.6 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.

5.7 Incubate the plate/tube at room temperature for ≥30 sec.

5.8 Carefully remove and discard the ethanol.

5.9 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.

5.10 Incubate the plate/tube at room temperature for ≥30 sec.

5.11 Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.

5.12 Allow the beads to dry at room temperature, sufficiently for all the ethanol to evaporate. **Caution: over-drying the beads may result in dramatic yield loss.**

5.13 Remove the plate/tube from the magnet.

6. Adapter Ligation Reaction Setup

- 6.1 To each well/tube containing the dried beads with A-tailed DNA, add:

Water (if required)	0 µl
Ligation Master Mix (Table 5A) [†]	45 µl
Adapter (concentration as required)	5 µl
Total volume per well/tube	50 µl

[†]For a safe stopping point, resuspend the beads in 20 µl 1X Ligation Buffer, cover the reaction and store at 4 °C for up to 24 hours. Resume the reaction by adding the rest of the reaction components, as outlined in Table 5B (p. 11).

- 6.2 Thoroughly resuspend the beads by pipetting up and down multiple times.
- 6.3 Incubate at 20 °C for 15 min. Proceed immediately to the next step.

7. Adapter Ligation Cleanup

Depending on your requirements and chosen workflow, either one or two post-ligation cleanups should be performed. Consult **Important Parameters: Post-ligation cleanup** (p. 5) and the **KAPA NGS Library Preparation Technical Guide** for more information.

First Post-Ligation Cleanup

- 7.1 To each 50 µl ligation reaction with beads, add:

PEG/NaCl SPRI® Solution	50 µl
Total volume per well/tube	100 µl

- 7.2 Thoroughly resuspend the beads by pipetting up and down multiple times.
- 7.3 Incubate the plate/tube at room temperature for 5 – 15 min to allow the DNA to bind to the beads.
- 7.4 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.
- 7.5 Carefully remove and discard the supernatant.
- 7.6 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.
- 7.7 Incubate the plate/tube at room temperature for ≥30 sec.
- 7.8 Carefully remove and discard the ethanol.
- 7.9 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.
- 7.10 Incubate the plate/tube at room temperature for ≥30 sec.
- 7.11 Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.

- 7.12 Allow the beads to dry at room temperature, sufficiently for all the ethanol to evaporate. **Caution: over-drying the beads may result in dramatic yield loss.**

- 7.13 Remove the plate/tube from the magnet.

- 7.14 Thoroughly resuspend the beads in an appropriate volume of elution buffer (10 mM Tris-HCl, pH 8.0). Consult **Important Parameters: Post-ligation cleanup** (p. 5). Recommended volumes are as follows:

- If proceeding to a second post-ligation cleanup, resuspend the beads in 50 µl.
- If proceeding to dual-SPRI® size selection, resuspend the beads in 100 µl, and omit the second post-ligation cleanup described below.
- For a safe stopping point, resuspend the beads in either 50 µl or 100 µl (as required), and store at 4 °C for up to 24 hours.

- 7.15 Incubate the plate/tube at room temperature for 2 min to allow the DNA to elute off the beads.

Second Post-Ligation Cleanup

- 7.16 To the 50 µl of cleaned-up, adapter-ligated DNA with beads, add:

PEG/NaCl SPRI® Solution	50 µl
Total volume	100 µl

- 7.17 Thoroughly resuspend the beads by pipetting up and down multiple times.
- 7.18 Incubate the plate/tube at room temperature for 5 – 15 min to allow the DNA to bind to the beads.
- 7.19 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.
- 7.20 Carefully remove and discard the supernatant.
- 7.21 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.
- 7.22 Incubate the plate/tube at room temperature for ≥30 sec.
- 7.23 Carefully remove and discard the ethanol.
- 7.24 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.
- 7.25 Incubate the plate/tube at room temperature for ≥30 sec.
- 7.26 Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.
- 7.27 Allow the beads to dry at room temperature, sufficiently for all the ethanol to evaporate. **Caution: over-drying the beads may result in dramatic yield loss.**

Second Post-Ligation Cleanup (continued)

- 7.28 Remove the plate/tube from the magnet.
- 7.29 Thoroughly resuspend the beads in 25 µl of elution buffer (10 mM Tris-HCl, pH 8.0). Please refer to the relevant notes under **Important Parameters**. Incubate the plate/tube at room temperature for 2 min to allow the DNA to elute off the beads.
- 7.30 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.
- 7.31 Transfer the clear supernatant to a new plate/tube and proceed with library amplification (Section 9), QC and/ or sequencing, as appropriate.

8. Dual-SPRI® Size Selection

- 8.1 To 100 µl of resuspended DNA with beads, add:

PEG/NaCl SPRI® Solution	60 µl
Total volume per well/tube	160 µl

- 8.2 Thoroughly resuspend the beads by pipetting up and down multiple times.
- 8.3 Incubate the plate/tube at room temperature for 5 – 15 min to allow library fragments larger than ~450 bp to bind to the beads.
- 8.4 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.
- 8.5 Carefully transfer 155 µl of the supernatant(s) containing library fragments smaller than ~450 bp to a new plate/tube. **It is critical to not transfer any beads with the supernatant.**
- 8.6 Discard the old plate/tube with the beads carrying library fragments larger than ~450 bp.
- 8.7 To each 155 µl supernatant in the new plate/tube, add:

Agencourt® AMPure® XP reagent	20 µl
Total volume per well/tube	175 µl

- 8.8 Thoroughly resuspend the beads by pipetting up and down multiple times.
- 8.9 Incubate the plate/tube at room temperature for 5 – 15 min to allow library fragments larger than ~250 bp to bind to the beads.
- 8.10 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.
- 8.11 Carefully remove and discard the supernatant.
- 8.12 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.
- 8.13 Incubate the plate/tube at room temperature for ≥30 sec.

- 8.14 Carefully remove and discard the ethanol.
- 8.15 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.
- 8.16 Incubate the plate/tube at room temperature for ≥30 sec.
- 8.17 Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.
- 8.18 Allow the beads to dry at room temperature, sufficiently for all the ethanol to evaporate. **Caution: over-drying the beads may result in dramatic yield loss.**
- 8.19 Remove the plate/tube from the magnet.
- 8.20 Thoroughly resuspend the beads in 25 µl of elution buffer (10 mM Tris-HCl, pH 8.0). Please refer to the relevant notes under **Important Parameters**. Incubate the plate/tube at room temperature for 2 min to allow the DNA to elute off the beads.
- 8.21 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.
- 8.22 Transfer the clear supernatant to a new plate/tube and proceed with library amplification (Section 9), QC and/or sequencing, as appropriate.

9. Library Amplification Reaction Setup

Please refer to **Important Parameters: Library amplification** and Section 7 of the **KAPA NGS Library Preparation Technical Guide** for more information on optimizing library amplification.

- 9.1 Assemble each library amplification reaction as follows:

Water (if required)	0 µl
Library DNA	20 µl
Library Amplification Master Mix (Table 6)	25 µl
PCR Primer Premix (5 µM each primer) [†]	5 µl
Total volume per well/tube	50 µl

[†]The recommended final concentration of each primer in the library amplification reaction is 500 nM. This may be increased to a final concentration of 2 µM of each primer if necessary.

- 9.2 Perform PCR with the thermocycling parameters outlined on the next page.
- 9.3 Store the plate/tube at 4 °C or -20 °C for up to 72 hours, or proceed directly to Section 10.

9. Library Amplification Reaction Setup (continued)

Step	Temp	Duration	Cycles
Initial denaturation	98 °C	45 sec	1
Denaturation	98 °C	15 sec	Minimum number required for optimal amplification (Table 2)
Annealing [†]	60 °C	30 sec	
Extension	72 °C	30 sec	
Final extension	72 °C	60 sec	1
Stop reaction	4 °C	Hold	Hold

[†]Optimization of the annealing temperature may be required for non-standard (i.e. other than Illumina® TruSeq™) adapter/primer combinations.

10. Library Amplification Cleanup

10.1 To each well/tube containing 50 µl library amplification reaction, add:

Agencourt® AMPure® XP reagent	50 µl
Total volume per well/tube	100 µl

10.2 Thoroughly resuspend the beads by pipetting up and down multiple times.

10.3 Incubate the plate/tube at room temperature for 5 – 15 min to allow the DNA to bind to the beads.

10.4 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.

10.5 Remove and discard the supernatant.

10.6 Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.

10.7 Incubate the plate/tube at room temperature for ≥30 sec.

10.8 Remove and discard the ethanol.

10.9. Keeping the plate/tube on the magnet, add 200 µl of 80% ethanol.

10.10 Incubate the plate/tube at room temperature for ≥30 sec.

10.11 Carefully remove and discard the ethanol. Try to remove all residual ethanol without disturbing the beads.

10.12 Allow the beads to dry at room temperature, sufficiently for all the ethanol to evaporate. **Caution: over-drying the beads may result in dramatic yield loss.**

10.13 Remove the plate/tube from the magnet.

10.14 Thoroughly resuspend the beads in an appropriate volume of elution buffer (10 mM Tris-HCl, pH 8.0). Please refer to the relevant notes under **Important Parameters**. Incubate the plate/tube at room temperature for 2 min to allow the DNA to elute off the beads.

10.15 Place the plate/tube on a magnet to capture the beads. Incubate until the liquid is clear.

10.16 Transfer the clear supernatant to a new plate/tube and proceed with library QC, target enrichment or sequencing, as appropriate.

Appendix 1:

Library construction guidelines for target enrichment with the Roche Nimblegen™ SeqCap EZ system

This Appendix provides guidelines for the preparation of libraries for targeted sequencing (capture) using the Roche Nimblegen™ SeqCap EZ system, and specific protocols for pre- and post-capture ligation-mediated (LM)-PCR. Please refer to the Roche Nimblegen™ SeqCap EZ SR User Guide (v3.0 and 4.0) for detailed instructions relating to hybridization of the amplified sample library to the SeqCap EZ library (capture), and washing and recovery of the captured DNA sample.

1. Library Construction

Prepare libraries from fragmented, human genomic DNA according to the protocol outlined on pp. 10 – 16 of this document. Specific recommendations are as follows:

- For optimal results with 2 x 100 bp paired-end sequencing, input DNA should be fragmented to an average size distribution of 180 – 220 bp, using a Covaris instrument. For 2 x 250 bp paired-end sequencing, a slightly larger insert size may be preferred. The size distribution of the fragmented DNA should be confirmed electrophoretically prior to library construction, using an PerkinElmer LabChip® GX, Agilent Bioanalyzer or TapeStation, or similar instrument.
- The input of fragmented, dsDNA into the end repair reaction should be in the range of 100 ng – 1 µg.
- The amount of adapter should be varied according to the amount of input DNA. Please refer to **Important Parameters: Adapter design and concentration** (p. 5) in this regard. When using Illumina® TruSeq™ adapter from an Illumina® TruSeq™ Library Preparation Kit v2:
 - For 1 µg input, dilute the adapter 1/5 in 10 mM Tris-HCl (pH 8.0) and use 5 µl per 50 µl ligation reaction.
 - For 100 ng input, dilute the adapter 1/50 in 10 mM Tris-HCl (pH 8.0) and use 5 µl per 50 µl ligation reaction.
- After adapter ligation, Roche typically recommends a dual-SPRI® size selection. This is achieved by performing the First Post-Ligation Cleanup (Steps 7.1 – 7.15), eluting the DNA in 100 µl elution buffer (10 mM Tris-HCl, pH 8.0), and proceeding to the Dual-SPRI® Size Selection (Section 8). The dual-SPRI® size selection is designed to eliminate adapter-ligated library fragments >450 bp and <250 bp. Please note that dual-SPRI® size selection can result in a significant loss of adapter-ligated material. If there is no need to eliminate large fragments (>450 bp) from the adapter-ligated library,

two consecutive 1X SPRI® cleanups (Steps 7.1 – 7.31) are recommended to remove unused adapter and adapter-dimer. This strategy results in a simpler workflow with less loss of material prior to pre-capture amplification, and should be considered in cases where input DNA is limiting.

- Irrespective of whether the "No Size Selection" (two consecutive 1X SPRI® cleanups) or Dual-SPRI® Size Selection strategy is followed, thoroughly resuspend the final, dried beads with cleaned/size-selected, adapter-ligated library in 25 µl elution buffer. Proceed directly to Section 2 below.

2. Pre-Capture LM-PCR

2.1 Assemble each pre-capture LM-PCR in a PCR tube or plate as follows:

2X KAPA HiFi HotStart ReadyMix	25 µl
5 µM TS-PCR Oligo 1	2.5 µl
5 µM TS-PCR Oligo 2	2.5 µl
Total volume per well/tube	30 µl

- 2.2 Place the plate/tube with beads resuspended in elution buffer (from Section 1) on a magnet to capture the beads. Incubate until the liquid is clear.
- 2.3 Transfer 20 µl of the clear supernatant (or PCR-grade water for negative controls) to the appropriate PCR tube or well with pre-capture LM-PCR mix. Mix thoroughly by pipetting up and down.
- 2.4 Perform the pre-capture LM-PCR with the following thermocycling parameters:

Step	Temp	Duration	Cycles
Initial denaturation	98 °C	45 sec	1
Denaturation	98 °C	15 sec	7 - 9†
Annealing	60 °C	30 sec	
Extension	72 °C	30 sec	
Final extension	72 °C	60 sec	1
Stop reaction	4 °C	Hold	Hold

†7 cycles are recommended for 1 µg input DNA. For lower input, or FFPE DNA, one or two additional cycles may be needed to produce sufficient material for capture.

2.5 Store the plate/tube at 4 °C for up to 72 hours or proceed directly to Section 3.

Please note:

- The sequences of the primers used for pre-capture LM-PCR may be found in the Roche Nimblegen™ SeqCap EZ SR User Guide v3.0.

2. Pre-Capture LM-PCR (notes, continued)

- The final concentration of each primer in the pre-capture LM-PCR is 250 nM. This should be increased to 500 nM when 1 µg input DNA is used for library construction, or if electrophoretic analysis of pre-capture LM-PCR products suggests that over-amplification is occurring. The best way to increase the final primer concentration is to increase the concentration of each primer stock from 5 to 10 µM, and use 2.5 µl of each primer.
- The forward and reverse primers may be supplied in a premix containing **each** of the primers at a final concentration of 2.5 µM (for a final concentration of 250 nM of each primer), or 5 µM (for a final concentration of 500 nM of each primer). If such a primer premix is used, the pre-capture LM-PCR mix will consist of 25 µl of KAPA HiFi HotStart ReadyMix (2X) and 5 µl of the primer premix.

3. Cleanup and QC of the Pre-Capture LM-PCR Amplified Sample Library

3.1 Amplified libraries can be purified using Agencourt® AMPure® XP reagent or the QIAGEN® QIAquick® PCR Purification Kit.

When using Agencourt® AMPure® XP reagent, follow the instructions in Section 10 (p. 16), with the following exceptions:

- Increase the volume of AMPure® XP reagent from 50 µl to 90 µl, for a 1.8X cleanup.
- Elute the cleaned library in 50 µl **PCR-grade water** instead of elution buffer.

When using the QIAGEN® QIAquick® PCR Purification Kit, follow the manufacturer's instructions, with the following exceptions:

- Elute the cleaned library in 50 µl **PCR-grade water** instead of Buffer EB.

3.2. Quantify the amplified library with the KAPA Library Quantification Kit, or with a NanoDrop™ spectrophotometer, and analyze 1 µl of the sample using an Agilent 2100 Bioanalyzer DNA 1000 Kit. For optimal results, the yield of amplified library should be >1 µg, and the library fragment size distribution should be 150 bp – 500 bp.

4. Hybridization of the Amplified Sample Library to the SeqCap EZ Library

4.1 Set up hybridization reactions according to the instructions in Chapter 5. **Hybridizing the Sample and SeqCap EZ Libraries** of the SeqCap EZ SR User Guide (v4.0).

5. Washing and Recovery of the Captured DNA Sample

5.1 Follow the instructions in Chapter 6. **Washing and Recovering Captured Multiplex DNA Sample** of the SeqCap EZ SR User Guide (v4.0) for binding of the captured DNA sample to streptavidin beads and washing of the bead-bound DNA.

6. Post-Capture LM-PCR

Two 50 µl post-capture LM-PCRs, each containing 20 µl of washed and resuspended streptavidin (capture) beads are performed for each library.

The composition of each post-capture LM-PCR is given in the table below. If convenient, double the amount of each component may be premixed (for a total of 60 µl of post-capture LM-PCR mix per library) and split equally between two PCR tubes or wells of a PCR plate.

Notes given at the end of Step 2.5 of the section on pre-capture LM-PCR also apply to post-capture LM-PCR.

6.1 Assemble **two** post-capture LM-PCRs for each captured library in PCR tubes or a PCR plate. Each reaction consists of:

2X KAPA HiFi HotStart ReadyMix	25 µl
5 µM TS-PCR Oligo 1	2.5 µl
5 µM TS-PCR Oligo 2	2.5 µl
Total volume per well/tube	30 µl

6.2 Vortex the streptavidin bead-bound captured Multiplex DNA Sample Library Pool to ensure that the solution is homogenous before proceeding to the next step.

6.3 Transfer 20 µl of streptavidin bead-bound capture DNA Sample Library (or PCR-grade water for negative controls) to each of the two wells/tubes containing post-capture LM-PCR mix. Mix thoroughly by pipetting up and down. Store the remaining bead-bound captured DNA Sample Library at -20 °C.

6.4 Perform the post-capture LM-PCR with the following thermocycling parameters:

Step	Temp	Duration	Cycles
Initial denaturation	98 °C	45 sec	1
Denaturation	98 °C	15 sec	14
Annealing	60 °C	30 sec	
Extension	72 °C	30 sec	
Final extension	72 °C	60 sec	1
Stop reaction	4 °C	Hold	Hold

6. Post-Capture LM-PCR (continued)

- 6.5 Store the plate/tube for up to 72 hours at 4 °C or proceed directly to Section 7.

7. Cleanup and QC of the Post-Capture LM-PCR Amplified Sample Library

- 7.1 Amplified libraries can be purified using Agencourt® AMPure® XP reagent or the QIAGEN® QIAquick® PCR Purification Kit.

When using Agencourt® AMPure® XP reagent, transfer the supernatants of each set of two PCRs from the streptavidin beads, and pool these into a single 1.5 ml microtube or prior to cleanup. Follow the instructions in Section 10 (p. 16), with the following exception:

- Increase the volume of AMPure® XP reagent to 216 µl, for a 1.8X cleanup.

Alternatively, the supernatants of the two corresponding PCRs may be cleaned up separately, and the DNA pooled after elution.

When using the QIAGEN® QIAquick® PCR Purification Kit, transfer the supernatants of each set of two PCRs from the streptavidin beads, and pool these into a single tube for purification using a single column. Follow the manufacturer's instructions and elute the final library in the appropriate volume of Buffer EB.

- 7.2 Quantify the final, amplified library with the KAPA Library Quantification Kit, or with a NanoDrop™ spectrophotometer, and analyze 1 µl of the sample using an Agilent 2100 Bioanalyzer DNA 1000 Kit. The yield of final, amplified library should be >500 ng, and the library fragment size distribution should be 150 bp – 500 bp.

8. Determine the Sample Enrichment using qPCR

- 8.1 Follow the instructions in **Chapter 8. Measuring Enrichment Using qPCR** of the **SeqCap EZ SR User Guide (v4.0)** for estimating the relative fold enrichment achieved for the amplified, captured DNA sample.

Restrictions and Liabilities

This technical data sheet is provided “as is” and Kapa Biosystems assumes no responsibility for any typographical, technical, or other inaccuracies. The document is subject to being changed, without notice, in future editions.

To the maximum extent permitted by applicable law, Kapa Biosystems disclaims all warranties, either express or implied, with regard to this technical data sheet and any information contained herein, including but not limited to the implied warranties of merchantability and fitness for a particular purpose. Kapa Biosystems shall not be liable for errors or for incidental or consequential damages in connection with the furnishing, use, or performance of this document or of any information contained herein.

This document might contain references to third party sources of information, hardware or software, products, or services and/or third party web sites (collectively the “Third-Party Information”). Kapa Biosystems does not control, and is not responsible for, any Third-Party Information. The inclusion of Third-Party Information in this document does not imply endorsement by Kapa Biosystems of the Third-Party Information or the third party in any way.

Kapa Biosystems is not responsible nor will be liable in any way for your use of any software or equipment that is not supplied by Kapa Biosystems in connection with your use of Kapa Biosystems products.

Kapa Biosystems does not in any way guarantee or represent that you will obtain satisfactory results from using Kapa Biosystems products as described herein. The only warranties provided to you are included in the Limited Warranty enclosed with this document. You assume all risk in connection with your use of Kapa Biosystems products.

Note to Purchaser: Limited Product Warranty

Any product that does not meet the performance standards stated in the product specification sheet will be replaced at no charge. This warranty limits our liability to the replacement of the product. No other warranties of any kind, express or implied, including without limitation, implied warranties of merchantability or fitness for a particular purpose, are provided by Kapa Biosystems. Kapa Biosystems shall have no liability for any direct, indirect, consequential or incidental damages arising out of the use, the results of use or the inability to use any product.

Note to Purchaser: Limited License

KAPA Library Preparation Kits are developed, designed and sold exclusively for research purposes and *in vitro* use. Neither the product, nor any individual component, has been tested for use in diagnostics or for drug development, nor is it suitable for administration to humans or animals. Please refer to the MSDS, which is available on request.

Certain applications of this product are covered by patents issued to parties other than Kapa Biosystems and applicable in certain countries. Purchase of this product does not include a license to perform any such applications. Users of this product may therefore be required to obtain a patent license depending upon the particular application and country in which the product is used.

Licensed under U.S. Patent nos. 5,338,671 and 5,587,287 and corresponding patents in other countries.

Agencourt®, AMPure® and SPRI® are registered trademarks of Beckman Coulter, Inc. Illumina® and TruSeq™ are trademarks or registered trademarks of Illumina, Inc. Nimblegen™ is a trademark of Roche. xGen™ and Lockdown™ are trademarks of Integrated DNA Technologies. Pippin Prep™ is a trademark of Sage Science. NanoDrop™ is a trademark of Thermo Fisher Scientific. LabChip® is a registered trademark of PerkinElmer, Inc. QIAGEN® and QIAquick® are registered trademarks of the QIAGEN group.

For support please contact: support@kapabiosystems.com

Boston, Massachusetts, United States
Int. Tel: +1 781 497 2933 Tel: +1 855 KAPA BIO
Fax: +1 781 497 2934
Email: sales@kapabiosystems.com

Cape Town, South Africa
Tel: +27 21 448 8200 Fax: +27 21 448 6503
Email: sales@kapabiosystems.com